

Product Information

MTSEA biotin-X

Catalog Number: M5036

Product Size: 10 mg

Application Scope: Thiol reactive dye


Parameters

Appearance: White solid soluble in DMSO or DMF

Molecular Formula: C₁₉H₃₄N₄O₅S₃

Molecular Weight: 494.68

Molecular Structural:

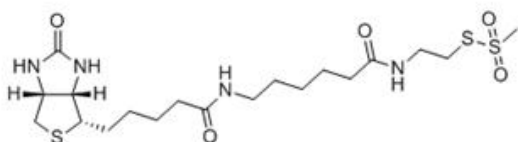

Storage

Store at desiccated -20°C and protect from light. When stored as directed, product is stable for at least 12 months.

Description

MTSEA biotin-X (MTSEA Biotincap; 2-((6-((biotinoyl) amino)hexanoyl)amino) ethylmethanethiosulfonate) can selectively label sulfhydryl groups on the cell surface under mild conditions.

Notes

1. There are quenching problems with fluorescent dyes. Please avoid light to slow down the fluorescence quenching.
2. For your safety and health, please wear lab coats and disposable gloves.

